


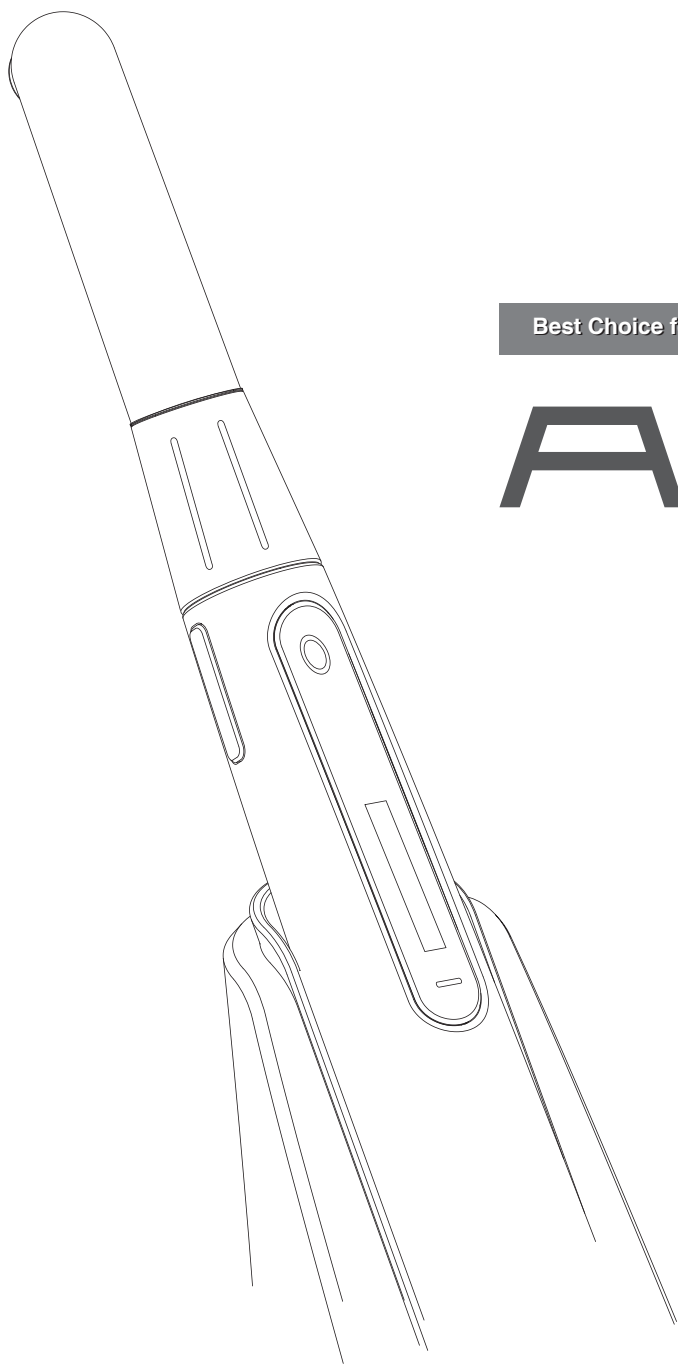

Best Choice for Best Treatment!

**A** Dr's Light **T**

Model : CL-AT24

User manual

Eng

CE

|                                              |    |
|----------------------------------------------|----|
| 1. Introduction                              | 04 |
| <hr/>                                        |    |
| - Features                                   |    |
| - Application                                |    |
| 2. Safety Information                        | 05 |
| <hr/>                                        |    |
| - Symbols                                    |    |
| - Notes and reference                        |    |
| - Installation and Connection                |    |
| - Operation                                  |    |
| - Maintenance                                |    |
| 3. List of Components                        | 09 |
| <hr/>                                        |    |
| 4. Parts and Functions                       | 10 |
| <hr/>                                        |    |
| 5. Installation and Connection               | 11 |
| <hr/>                                        |    |
| 6. Direction for Use                         | 13 |
| <hr/>                                        |    |
| 7. Selection of Curing Modes                 | 14 |
| <hr/>                                        |    |
| 8. Maintenance                               | 15 |
| <hr/>                                        |    |
| - Cleaning and sterilization                 |    |
| - Maintenance and replacement of the Battery |    |
| 9. Trouble Shooting                          | 16 |
| <hr/>                                        |    |
| 10. Specification                            | 17 |
| <hr/>                                        |    |
| 11. Accessories                              | 18 |
| <hr/>                                        |    |
| 12. Warranty                                 | 18 |
| <hr/>                                        |    |
| 13. Warranty card                            | 19 |
| <hr/>                                        |    |

# 1 Introduction

Thank you for purchasing Dr's Light AT.

The instructions described below have been designed to provide necessary information for successful operation, daily maintenance and safety for use of this unit.

## Features

---

- Variety of curing modes
- Easy-to read LCD screen
- Enhanced capacity battery
- Lithium-ion battery
- Interchangeable lithium-ion battery
- Chargeable battery alone
- Automatic sleep mode


## Application

---

- Polymerization of light-cure catalyst containing dental materials with a wavelength spectrum of 440 ~ 490 nm.

## Overview

---

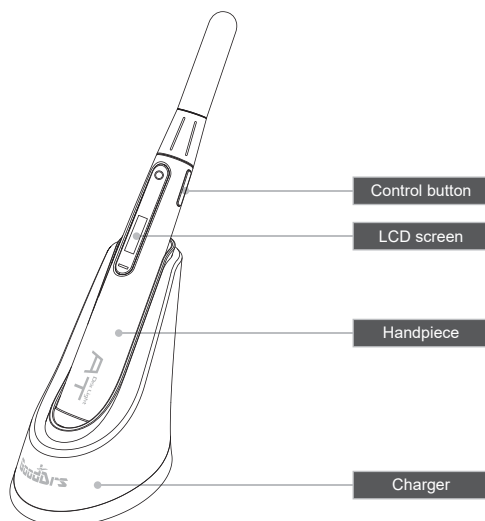

## 2 Safety information

It is important to understand the following information to use this product safely.

### Symbols

---



Indicates a potentially hazardous situation which, if not avoided, may result in serious injury or death.



Indicates a potentially hazardous situation which, if not avoided, may result in minor/moderate injury or in damage of property

### Notes and References

---

Notes and References used in this manual are shown below.



Indicates that potential malfunctions may occur if the description is ignored.



Supplementary information to facilitate better operation.

## 2 Safety information

### Installation and Connection



#### Warning

- The adaptor should be disconnected from the charger before connecting the power cable to prevent an electric shock.
- Keep away from sources of ignition to avoid fire or explosion.
- Protect the handpiece, charger, and battery from any liquid to avoid electric shock or functional error.
- Place the charger on a flat and stable surface. Do not drop handpiece or charger in order to avoid internal defect that may cause overheating or fire.
- Do not place the electrical cord near a gas burner or sources of heat.  
Should the cord catch on fire, replace with a new one to avoid short-circuit.



#### Caution

- Do not use or leave the unit in areas of high temperature or exposed to direct sunlight; internal circuit defect may cause overheating or fire.
- Securely connect the dust-free power plug to the outlet.  
Use this device in well ventilated areas. Overheating may cause fire.
- Do not touch the terminal areas of charger or battery with fingers or conductive objects, such as metal wire or safety-pin. It may cause electric shock or fire.
- Do not use any charger other than supplied. It may result in electric shock or fire.

## 2 Safety information

### Operation



#### Warning

- Do not plug the power cord with wet hands. This may result in an electric shock.
- In case of fluid leak from battery or partial discoloration/deformation of the handpiece, immediately discontinue use.
- Place the battery in correct direction when charging. Wrong directed placement may cause malfunction, explosion, or fire.
- In case of partial discoloration of the charger, immediately disconnect the power plug and repair the charger.
- In case battery fluid contacts the eye, immediately flush with copious amount of clean water and seek medical attention.
- Do not use this product on patients with a history of cataract, retinal damage or allergic reaction to light.
- Do not use this product on a person with a pacemaker as it may cause malfunction of the pacemaker.
- Avoid exposing the light on intraoral soft tissue as it may affect the soft tissue.
- Do not put the battery on the charger upside down. This may cause an explosion or fire.



#### Caution

- Disconnect the power cord by grasping the power plug to avoid injuries or burns.
- Avoid damage or use of unnecessary force on the power cord. This may cause an electric shock or fire.
- Should battery fluid contact the skin or clothes, immediately wash it off with copious amounts of water to avoid skin damage.
- Take the battery out of the handpiece if it won't be used for more than 1 week.
- Connect the power plug only when battery is being charged. Insulation degradation may cause electric shock or fire.
- Do not use this product for any purposes other than dental use.

## 2 Safety information

### Maintenance

---



#### Warning

- Use only the battery supplied to avoid fluid leak or explosion.
- Never disassemble or overhaul the product. Unusual performance may cause injuries, electric shock, or fire.
- Ensure that charger and battery are always separated from the handpiece prior to cleaning.




#### Caution


- If the LCD screen malfunctions, remove the battery from the handpiece and reinsert it.

### 3 List of components


Please check the contents described below when removing the contents from the original packaging.

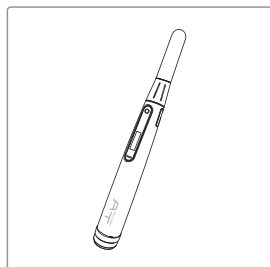


Handpiece

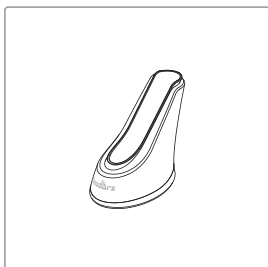


Charger

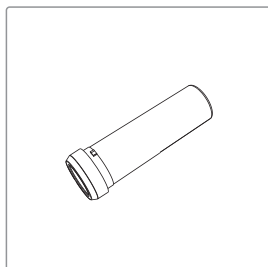


Battery

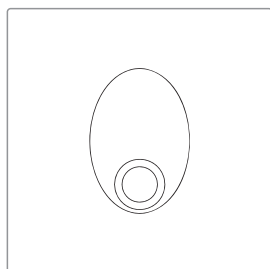


Anti-glare shield

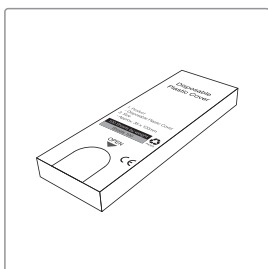

Disposable cover  
(Approx. 100pcs)

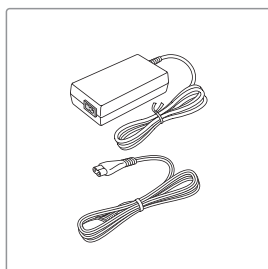


Adapter & Power cable

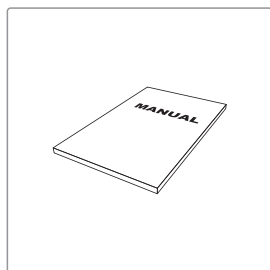

User manual

## 4 Parts and functions

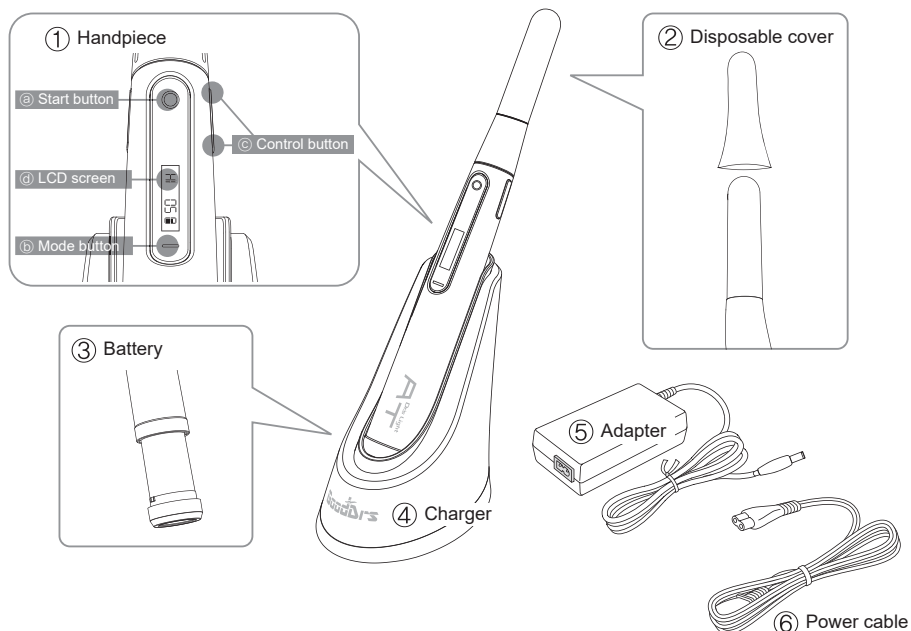

### ① Handpiece

- a Start button : Operation and Turn On/Off
- b Mode button : Curing mode select
- c Control button : Curing time select
- d LCD screen : Display device state

### ② Disposable cover (Approx.100pcs)

### ③ Battery

### ④ Charger : Charge battery. Display charging process

### ⑤ Adapter


### ⑥ Power cable

## 5 Installation and Connection

Refer to the information supplied on safety precautions during installation and connection.

### 1. Connection between Adapter and Charger


- 1) Securely insert the power cable with 100~240V power outlet.
- 2) Connect power cable with adapter.
- 3) Connect adapter with charger

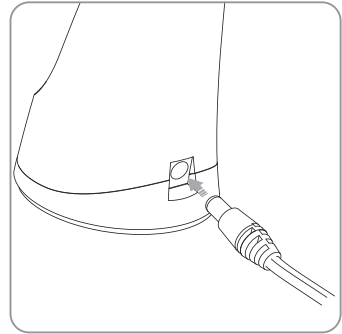

- Use only the original adapter supplied.
- When the adapter is connected with charger, red and blue light will blink.
- If red and blue light is not blinking when adapter is connected, it is a charger malfunction.

### 2. Inserting Battery


- ① Carefully insert the battery into the handpiece as directed by the arrow.
  - ② The battery turns with the arrow, as seen in image.
- \*Avoid forcing the Battery.

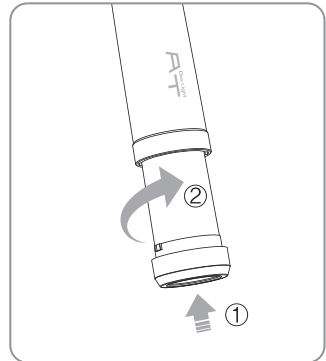

## 5 Installation and Connection

### 3. Charging

- Insert the handpiece into the vertical slot at the front of the charger.
- Red light on : charging in progress
- Blue light on : charging process is complete

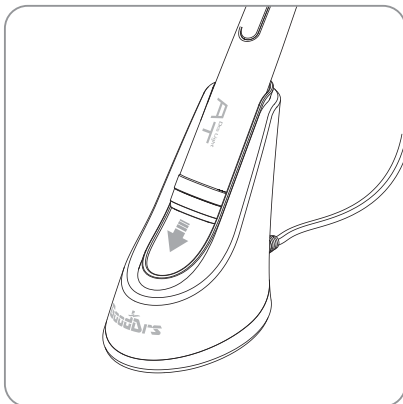

Note

- If the remaining battery is extremely low, red light may flicker while recharging, which is not a malfunction.



Reference

- It will take approximately 3~4 hours before the battery will be fully charged.


## 6 Direction for Use

### 1. Activating Handpiece



- 1) Press any button, then LCD screen will be on standby mode and time display used last.
- 2) Select curing mode and time setting with mode button and control button.



Initial Screen

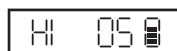

Standby Screen  
(last used mode & time)

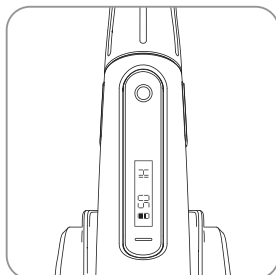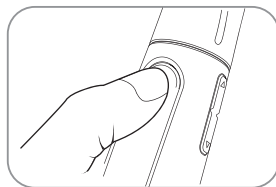

Reference

- When Battery inserted, the unit is automatically turned on.

### 2. Start of exposure

Press Start button to start exposure with the selected mode.



Reference

- Press any button if you wish to end the curing process.
- Alarm during operating informs time passage.

### 3. Battery icon



This icon indicates the power remaining in the battery. Recharge or replace the battery when "LOW BATTERY" sign appears on the display screen.



Reference

- Sleep Mode : If unit is left unused for a few minutes, all functions are automatically transferred to sleep mode with a short beep. Hold down start button for approximately 3-5 seconds to shut off light. If you want to turn on, press any button.
- Set-up memory : The handpiece will record the setting used last.


## 7 Selection of Curing Mode

Curing mode can be set by pressing the mode button.  
Various curing programs are available in this unit.

### 1. High Mode

Time setting: 5, 10 sec.


Quick polymerizing with high intensity

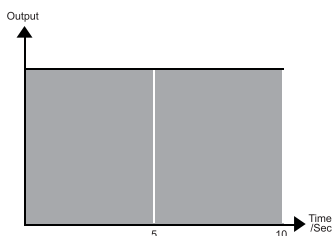

### 2. Low Mode

Time setting: 10, 20, 30 sec.


Polymerizing with low intensity

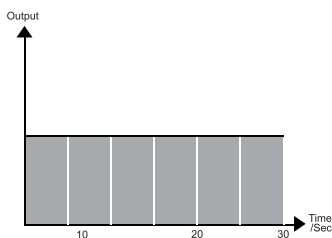

### 3. Soft Start Mode

Time setting: 20 sec. (fixed)


Special mode to avoid polymerization shrinkage

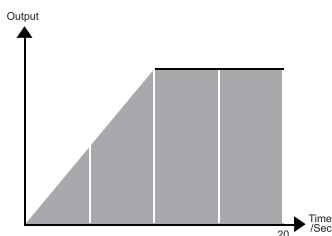

### 4. Pulse Mode

Time Setting: 10, 20 sec.

Special mode to avoid polymerization shrinkage

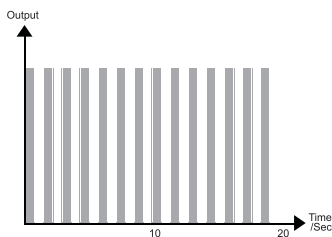

### **Cleaning and disinfection**

---

- Clean the handpiece and charger gently with clean, soft, lint-free cloth.
- Avoid using solvents which contain ethanol or acetone.
- Do not use sharp or pointed tools to clean as they may damage or scratch the plastic.
- Be careful not to spill any liquid into the device.

### **Maintenance and replacement of the battery**

---

- Exposure to high room temperature may shorten the battery life and affect device performance. Use the device in normal temperature.
- It is recommended for rechargeable battery to be charged immediately before use to maximize efficiency.
- Remove the battery from the handpiece, if not being used for a long time.
- If the running time declines with a full-charged battery, replace the battery before use.

## 9 Trouble Shooting

If your curing light unit is not operating normally, please refer to the following points.

| Trouble                      | Cause and Solution                                                                                                                                                                                                                            |
|------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Control button does not work | If the control button is pressed with excess force, it will not work. In such case press the button gently.                                                                                                                                   |
| Display screen does not work | If the battery is completely discharged or connected incorrectly, display will not appear. Recharge the battery or reconnect the battery to reset the device.                                                                                 |
| Power off while in operation | <p>If the start button is pressed for a moment, power will be turned off. Press any button again. Check the battery connection again.</p> <p>If the battery is completely discharged, device can be turned off.<br/>Recharge the battery.</p> |
| Not able to charge           | If the power cord or etc. is not correctly connected, the battery can not be charged. If the battery or handpiece is not properly placed on the charger, the battery can not be charged.                                                      |

\* If the trouble is not settled after the above solutions, please contact local dealer or Good Doctors Co., Ltd.

## 10 Specification

|           |                               |                          |
|-----------|-------------------------------|--------------------------|
| Handpiece | Operation voltage             | 3.7 V                    |
|           | Wavelength range              | 440 ~ 490 nm             |
|           | Dimensions (included Battery) | Approx. 232 x 26 x 26 mm |
|           | Weight (included Battery)     | Approx. 130 g            |

|         |                   |                           |
|---------|-------------------|---------------------------|
| Charger | Operation voltage | 5 V                       |
|         | Dimensions        | Approx. 110 x 105 x 68 mm |
|         | Weight            | Approx. 100 g             |

|         |                 |                         |
|---------|-----------------|-------------------------|
| Battery | Type            | Lithium Ion Battery     |
|         | Nominal Voltage | DC 3.7 V                |
|         | Dimensions      | Approx. 75 x 25 x 25 mm |
|         | Weight          | Approx. 52 g            |

|                  |            |                           |
|------------------|------------|---------------------------|
| Disposable cover | Dimensions | Approx. 32 x 100 mm       |
|                  | Quantity   | Approx. 100 pcs by weight |


|         |                |                |
|---------|----------------|----------------|
| Adapter | Input voltage  | AC 100 ~ 240 V |
|         | Output voltage | DC 5V          |

### ■ Operation, transport, and storage condition for the main unit


|                      | Operation      | Storage | Transport |
|----------------------|----------------|---------|-----------|
| Temperature          | 10~40C         | 0~40C   | -10C~50C  |
| Humidity             | 10~80%         | 0~80%   | 0~80%     |
| Atmospheric pressure | 800 ~ 1060 hPa |         |           |

\* This specification is subject to change without prior notice


## 11 Accessories

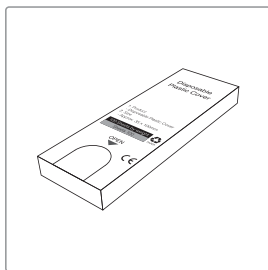

Disposable cover  
(Approx.100pcs)

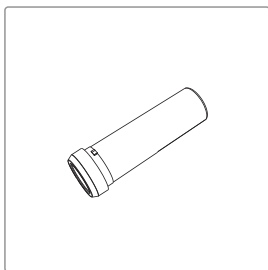

Battery

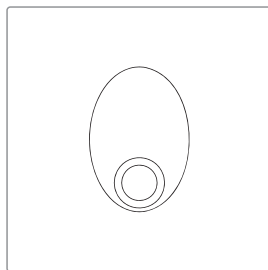

Anti-glare shield



## 12 Warranty

This product is shipped after strict inspection. Should the product have any problem under normal usage within one year from the date of purchase, repair will be free of charge.

However, repair of consumables (disposable cover, battery, etc.) are not covered by the warranty.

For more information, refer to the warranty card attached.

### Serial number check.

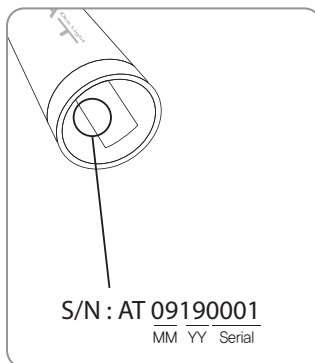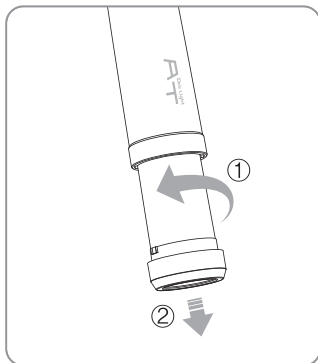

### Good Doctors Co., Ltd.

#208, B-dong, 283 Bupyeong-daero, Bupyeong-gu  
(Woolim Lions Valley, Cheongchoen-dong), Incheon,  
21315, Republic of Korea  
Tel : + 82 32 424 6325 Fax : + 82 32 424 6326  
E-mail : info@gooddrs.com Website : www.gooddrs.com

### EC REP Good Doctors Germany GmbH

Gerhard-Domagk-Str. 2, 53121 Bonn, Germany  
TEL: +49 (0) 228/53 44 14 65 FAX: +49 (0) 228/53 44 14 66  
www.gooddrs.de



We guarantee the quality of the product as follows

**Name of product :**

**Name of model :**

**Customer information**

- Office name :
- Office phone :

**Distributor**

- Dealer name :

**Assurance period**

- Date of purchase : from , , 20
- Product is for 1 year except accessories.

**< Assurance details >**

1. If there is any defect of manufacturer or product has any fault itself within the assurance period, it will be replaced with the same product.
2. Below cases are not effective for the free refund.
  - Malfunction and damage according to user's carelessness.
  - Malfunction and damage according to natural disaster such as fire, earthquake and flood etc.
  - Malfunction and damage according to user's intention or accident.
  - Malfunction and damage according to using incompatible product.
3. The manufacturer does not accept any liability resulting misuse.



**Good Doctors Co., Ltd.**

#208, B-dong, 283 Bupyeong-daero, Bupyeong-gu (Woolim Lions Valley, Cheongchoen-dong) Incheon 21315 REPUBLIC OF KOREA

Tel: +82 32 424 6325 Fax: +82 32 424 6326